


# MC ETD2 Wrist Disarticulation

## Prosthetist Manual

### Introduction

The Wrist Disarticulation option provides an attachment plate that can be laminated directly to the outer socket of the prosthesis. This minimizes the length of the prosthesis in cases where much of the patient's remnant limb has been saved. While it does allow a slight amount of rotational adjustment, it sacrifices passive rotation and interchangeability of terminal devices for the patient.

### Components

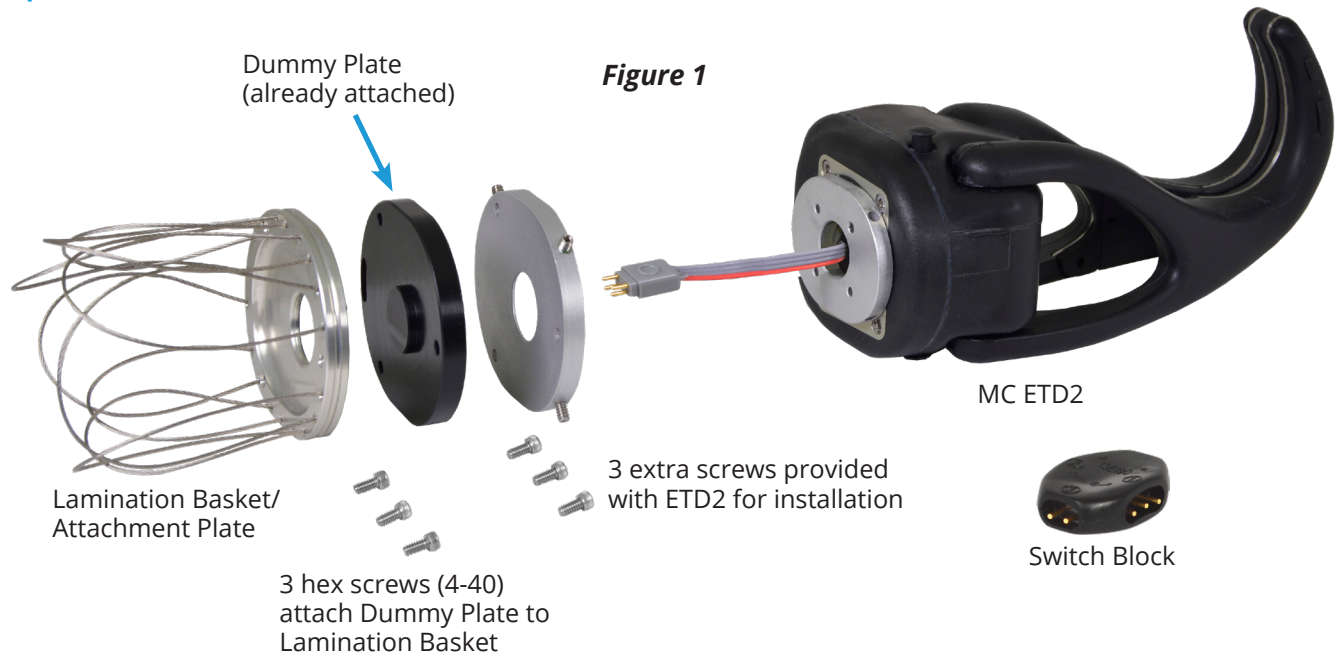

### Instructions

1. Using foam, plaster of Paris, wax or other preferred material, shape the distal portion of the forearm. Ensure that the inner socket can be removed from the outer socket to allow installation of electrodes, preamps, wires, and internal battery.
2. Establish the mounting surface for the Lamination Basket. This must be at the correct length and wrist angle.
3. If necessary, use a parting agent (such as a PVA bag) between the forearm material and the Lamination Basket. This will allow easy removal of material, and removal of the inner socket from the outer forearm.
4. Temporarily attach the Lamination Basket, making sure the "TIPS" marking on the Lamination Dummy Plate (already installed) is in the direction of the MC ETD2 fingertips (Figure 2). Otherwise the lamination basket will be off by 180 degrees.
5. In the lay up, make sure at least 2 stockinettes (4 layers) are firmly tied off in the Tie-Off Grove.
6. Remove and add parting agent to each of the 3 hex screws holding the lamination dummy in place.
7. Reinstall the Lamination Dummy Plate and screws.
8. Laminate.
9. Trim the excess lamination away and remove the screws and the lamination dummy. Use care in removing the last piece of plastic to prevent damage to the attachment plate.

10. With the inner socket removed, feed the wire of the MC ETD2 through the center slot in the attachment plate.
11. Place a small amount of medium strength threadlocker (Loctite 242, 243 or equivalent) on the 3 hex screws packaged with the MC ETD2. Use these screws to attach the MC ETD2 to the attachment plate. Note that the attachment holes are slotted to allow some rotational adjustment.
12. With a 3/32 in. ball end hex key, firmly hand-tighten the screws.



**Figure 2**

## Other Information

Please see the enclosed MC ETD2 Prosthetist Manual and User Guide for important Instructions, Cautions, Specifications, Warranty information, etc.

## Serious Incidents

In the unlikely event a serious incident occurs in relation to the use of the device, users should seek immediate medical help and contact their prosthetist, local competent authority and Fillauer at the earliest possible convenience. Clinicians should at any time contact their local Fillauer representative and local competent authority immediately in the event of any device failure.

## Declaration of Conformity

The product herewith complies with the Medical Device Directive 93/42/EEC guidelines and is registered with the United States Food and Drug Administration (Registration No. 1723997).



## Customer Support

### Americas, Oceania, Japan

#### Fillauer Motion Control

115 N. Wright Brothers Drive  
Salt Lake City, UT 84116  
801-326-3434  
motioninfo@fillauer.com

### Europe, Africa, Asia

#### Fillauer Europe

Kung Hans väg 2  
192 68 Sollentuna, Sweden  
+46 (0)8 505 332 00  
support@fillauer.com

*Fillauer*®

www.fillauer.com



#### Fillauer

2710 Amnicola Highway  
Chattanooga, TN 37406  
423.624.0946



#### Fillauer Europe

Kung Hans väg 2  
192 68 Sollentuna, Sweden  
+46 (0)8 505 332 00



© Fillauer Motion Control, 1910106 Rev D 02-18-25